



**Surflo®**

**Winged infusion sets**




**Smooth  
vein access  
Better comfort  
for the patient**

**TERUMO®**

We keep life flowing


## Unsurpassed penetration performance

Surflo® Winged infusion sets incorporate the best of TERUMO®'s needle technology. The precisely honed edges and triple bevel design facilitate the smoothest possible penetration with minimum tissue trauma. The lubricant treatment on the cannula guarantees trouble free insertion into the vein.

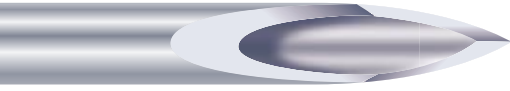

## Ultra thin-walled for higher flow rates

The ultra thin-walled cannula allows for higher flow rates and therefore is ideal for rapid infusions and dosages with high viscosity.

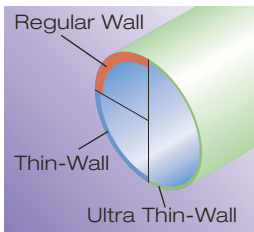

## Flexible, memory-less tubing and transparent luer adaptor

The material used for the tubing is highly flexible and virtually memory-less and therefore eliminates the risk of kinking and disruption of the infusion. For added convenience, both the tubing and the luer adaptor are translucent for easy observation of the contents.




## 27G dedicated to neonatology



## Large colour coded wings

Surflo® winged infusion sets are equipped with large, interlocking wings which ensure perfect stability and grip during venipuncture. To conform to the skin surface and to put less pressure on the blood vessel, the wings are slightly angled in relation to the cannula. The wings are colour-coded for easy gauge identification.

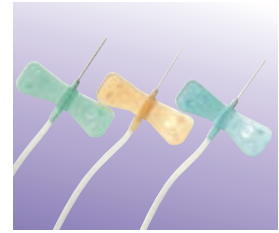

TERUMO® references	EN/ISO colour code	External diameter of the cannula		Cannula length mm	Dead space volume ml	Minimum flow rate ml/min	Tubing length mm
		Gauge	mm				
SV-18BL	Pink	18	1.2	19	0.44	47	300
SV-19BL	Ivory	19	1.1	19	0.43	32	300
SV-21BL	Green	21	0.8	19	0.42	13.5	300
SV-22BL	Black	22	0.7	19	0.42	8.2	300
SV-23BL	Blue	23	0.6	19	0.41	2.9	300
SV-25BL	Orange	25	0.5	19	0.41	1.6	300
SV*27EL	Grey	27	0.4	13	0.27	1.5	200



Every day, as healthcare professionals, you are working to ensure life and to contribute to the well-being of your patients. Every day, Terumo Europe provides you with the clinical systems, means and solutions to keep life flowing.



**TERUMO®**

We keep life flowing

### TERUMO CORPORATION

44-1, 2-chome  
Hatagay  
Shibuya-ku  
Tokyo 151-0072  
JAPAN  
Tel. +81/3 3374.8111  
Fax +81/3 3374.8196

### TERUMO EUROPE N.V.

Researchpark Zone 2 Haasrode  
Interleuvenlaan 40  
B-3001 Leuven  
BELGIUM  
Tel. +32/16.38.12.11  
Fax +32/16.40.02.49

### TERUMO EUROPE N.V. BENELUX SALES DIVISION

Researchpark Zone 2 Haasrode  
Interleuvenlaan 40  
B-3001 Leuven  
BELGIUM  
Tel. +32/16.39.25.80  
Fax +32/16.39.25.99  
THE NETHERLANDS  
Tel. 0800/022.03.96  
Fax 0800/022.04.14

### LABORATOIRES TERUMO

FRANCE S.A.  
1, Parc Ariane  
Bâtiment Uranus  
Rue Hélène Boucher  
78284 Guyancourt Cedex  
FRANCE  
Tel. +33/1 30.96.13.00  
Fax +33/1 30.43.60.85

### TERUMO DEUTSCHLAND GmbH

Hauptstrasse 87  
65760 Eschborn  
GERMANY  
Tel. +49/61.96.80.230  
Fax +49/61.96.80.23.200

### TERUMO EUROPE ESPAÑA SL

Avenida de Burgos 16-D  
Edificio EUROMOR  
Madrid 28036  
ESPAÑA  
Tel. +34/91 383.90.69  
+34/91 383.91.61  
Fax +34/91 767.04.60

### TERUMO ITALIA SRL

Via Simone Martini 143/145  
I-00142 Roma  
ITALIA  
Tel. +39/06.519.61.420  
Fax +39/06.503.04.07

### TERUMO SWEDEN AB

Nya Varvet  
Byggnad 90  
SE-426 71 Västra Frölunda  
SWEDEN  
Tel. +46/31.748.58.80  
Fax +46/31.748.58.90

### TERUMO DENMARK Filial af TERUMO SWEDEN AB

Tel. +45/7020.93.80  
Fax +45/7020.94.80

### TERUMO UK Ltd.

Tamesis  
The Causeway  
Egham  
Surrey  
TW20 9AW  
UNITED KINGDOM  
Tel. +44/1784.476.200  
Fax +44/1784.476.201

®: Registered Trademark  
Published by Terumo Europe N.V.

[www.terumo.com](http://www.terumo.com)  
[www.terumo-europe.com](http://www.terumo-europe.com)